

CARD16 Rabbit pAb

Catalog No.	A16624	Category	Polyclonal Antibodies
Applications	WB	Observed MW	28kDa
Cross-reactivity	Mouse	Calculated MW	10kDa/22kDa

Immunogen Information

Immunogen	A synthetic peptide corresponding to a sequence within amino acids 50-150 of human CARD16 (NP_001017534.1).
Gene ID	114769
Swiss prot	Q5EG05
Synonyms	CARD16;COP;COP1;PSEUDO-ICE


Product information

Source	Rabbit
Isotype	IgG
Purification method	Affinity purification
Storage	Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide, 50% glycerol, pH 7.3.

Background

Recommended Dilutions

WB 1:500 - 1:2000

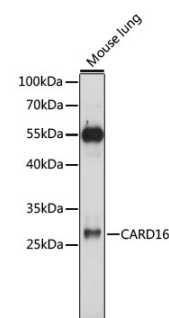

Western blot analysis of extracts of mouse lung, using CARD16 antibody (A16624) at 1:1000 dilution.
Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25ug per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 90s.