

NUAK1 Rabbit pAb

Catalog No.	A17333	Category	Polyclonal Antibodies
Applications	WB	Observed MW	74kDa
Cross-reactivity	Human, Rat	Calculated MW	23kDa/74kDa

Immunogen Information

Immunogen	Recombinant fusion protein containing a sequence corresponding to amino acids 362-661 of human NUAK1 (NP_055655.1).
Gene ID	9891
Swiss prot	O60285
Synonyms	NUAK1;ARK5


Product information

Source	Rabbit
Isotype	IgG
Purification method	Affinity purification
Storage	Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide, 50% glycerol, pH 7.3.

Background

Recommended Dilutions

WB 1:200 - 1:2000

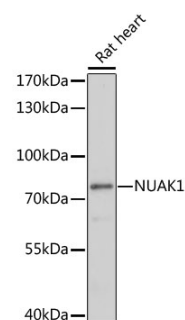

Western blot analysis of extracts of rat heart, using NUAK1 antibody (A17333) at 1:1000 dilution.
Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25ug per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 30s.